


1 Electrodes



channel 1 FCR · channel 2 ECR · one reference, on the olecranon

2 Connect

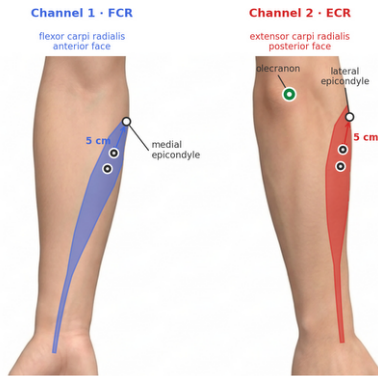

with the board switched on · in «Test identifier», a code, never a name

3 Record

Start recording

Warm up first

10

Two or three easy contractions of each muscle. The first maximal effort of a session is never the strongest one.

warm-up: 10 s, two or three easy contractions

4 Calibration



Each effort: 3 s warning → one brief, explosive jerk within the 1.5 s count → 2 s rest



When it says FCR

a maximal jerk of wrist flexion, the fist clenched with all your strength




When it says ECR

a maximal jerk of wrist extension, hand open, fingers stretched out

3 × FCR · 3 × ECR · one explosive jerk in the 1.5 s of each cue

5 The task

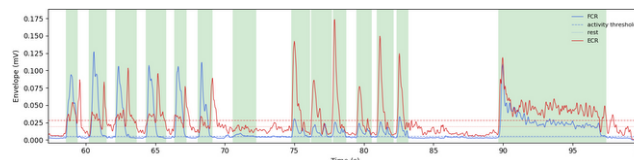

Stop recording

6 flexions · 6 extensions · 1 grip of 8 s squeezing the ball · 2 s still in between

6 Analyse

Select fragments...

Keep it

Drop it

Split it

Use these fragments

Generate PDF report

the yellow line says the step; last, the PDF report

If it does not connect: if the board does not answer, wait for the warning and press «Retry»; if Windows gives an error opening the port, switch the board off and on and press «Start recording» again.

In the analysis: check that every contraction has its row; the very brief ones may come out dotted or unmarked, and are added by hand (clicking the dotted one or with «Add fragment»).

Questions? The «?» on each box of the application explains it.