


1 Electrodes



channel 1 FCR · channel 2 ECR · reference on the ulna

2 Connect

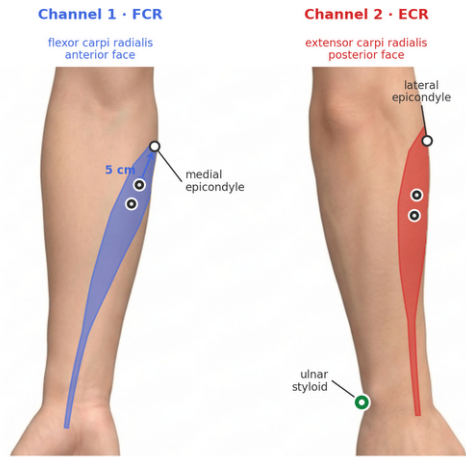

with the board switched on

3 Record

Start recording

Warm up first

10

Two or three easy contractions of each muscle. The first maximal effort of a session is never the strongest one.

warm-up: 10 s, two or three easy contractions

4 Calibration




Each effort: 3 s warning → one brief, explosive jerk within the 1.5 s count → 2 s rest



3 × FCR · 3 × ECR · one explosive jerk in the 1.5 s of each cue

5 The task

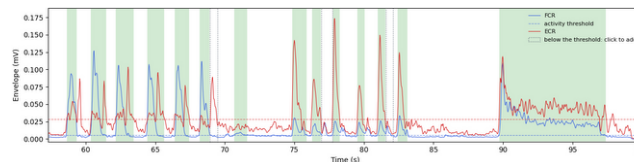

Stop recording

6 flexions · 6 extensions · 1 grip of 5 s · 2 s still in between

6 Analyse

Select fragments...

Keep it

Drop it

Split it

Use these fragments

Generate PDF report

the yellow line says the step; last, the PDF report