


1 Electrodes



channel 1 FCR · channel 2 ECR · one reference, on the olecranon

2 Connect

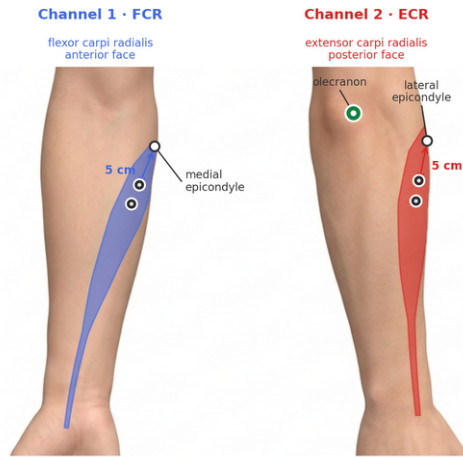

with the board switched on · in «Test identifier», a code, never a name

3 Record

Start recording

Warm up first

10

Two or three easy contractions of each muscle. The first maximal effort of a session is never the strongest one.


warm-up: 10 s, two or three easy contractions

4 Calibration



3 × FCR · 3 × ECR · one explosive jerk in the 1.5 s of each cue

5 The task

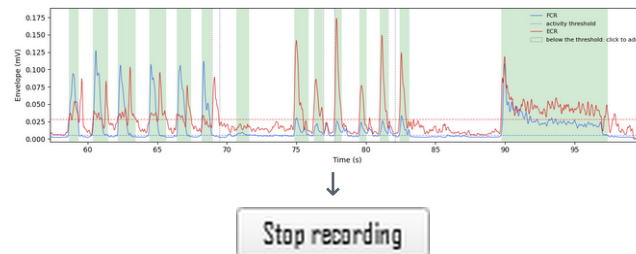

Stop recording

6 flexions · 6 extensions · 1 grip of 8 s · 2 s still in between

6 Analyse

Select fragments...

Keep it

Drop it

Split it

Use these fragments

Generate PDF report

the yellow line says the step; last, the PDF report