


1 Electrodes



channel 1 FCR · channel 2 ECR · one reference, on the olecranon

2 Connect

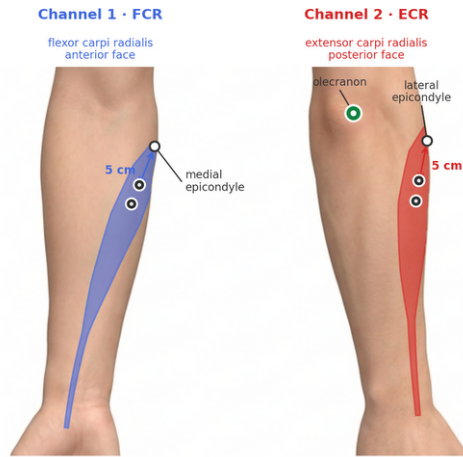

with the board switched on · in «Test identifier», a code, never a name

3 Record

Start recording

Warm up first

10

Two or three easy contractions of each muscle. The first maximal effort of a session is never the strongest one.


warm-up: 10 s, two or three easy contractions

4 Calibration




3 × FCR · 3 × ECR · one explosive jerk in the 1.5 s of each cue

5 The task

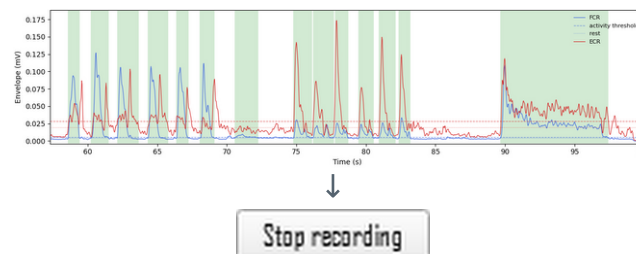

6 flexions · 6 extensions · 1 grip of 8 s squeezing the ball · 2 s still in between

6 Analyse

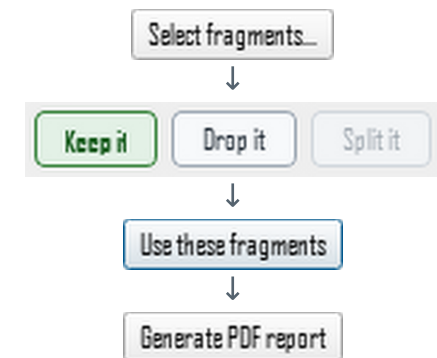

the yellow line says the step; last, the PDF report